

BP37™




BP37™ Infant Temperature Maintenance Kit is designed to provide a consistent and uniform surface temperature of 37°C for the infant to lie on moments after switching it on.

Simple to use, the single use Thermal Support Matrix is easily pushed into the Clamp (TSMC) and the Controller is turned on. Pre-set to 37°C, the Matrix will warm up in seconds to allow the baby to be placed directly onto a warm surface to help them maintain their temperature. Low tack adhesive tabs are used to secure the Matrix to a mattress, ensuring that it will not slide when the baby is placed on it during transport.

The lightweight BP37™ Controller clearly displays actual surface temperature, ambient temperature, battery percentage and time since activation. It is supplied with a fixing kit to enable the Controller to be attached to AHT Baby Pod® range, optional clamps allow the Controller to be mounted onto a Medi-Rail or a standard Pole. Running from internal batteries for up to 6 hours, it can be fully recharged in approximately 3 hours, allowing flexibility to be used away from mains electricity for prolonged periods. When thermal support is not needed, the Controller can be switched to standby mode to preserve the battery.

Technical Details	
Dimensions (LxWxH)	125 x 115 x 95mm
Weight (Controller)	610g (1.3lbs)
Controller Materials	Body - PA66
	Screen - Stalinate
	TSMC Cable - PVC
	TSMC Clamp - PA66
Matrix Materials	Top surface - Non-woven Spun bond PP
	Adhesive Pads - Acrylic solvent on a PE carrier
	Backing - Polyester
	Printed Circuit - Aluminium
Safety Compliance	CE Mark Class I
	FDA Device Class 2
Battery Usage Time	Up to 6 hours
Battery Recharge Time (0 to 100%)	Up to 3 hours

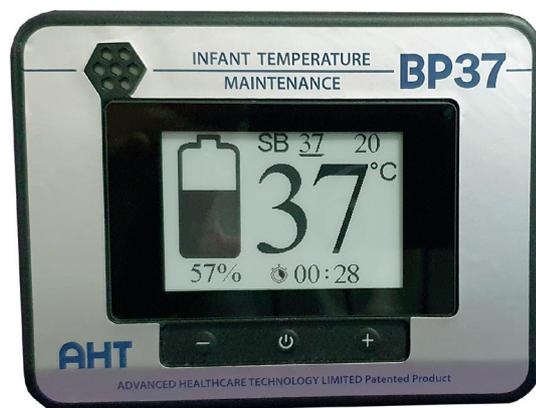

For more information on our products, visit our website at: www.babypod.com
or E-mail us at: sales@babypod.com

Ordering Information

Product Description	
Part Number	Product Description
BP37-001	BP37™ Infant Temperature Maintenance Kit Starter Pack Complete with Controller, Battery Charger & Power Cord, TSMC, Fitting Kit and single use Matrix (Pack of 10)
BP37-001-LS	BP37™ Infant Temperature Maintenance Kit Starter Pack - LS Variant Complete with Controller, Battery Charger & Power Cord, TSMC, Rail Clamp, Cable Guide and single use Matrix (Pack of 10)
BP37-002	BP37™ Controller Complete with Controller, Battery Charger & Power Cord, TSMC & Fitting Kit
BP37-002-LS	BP37™ Controller - LS Variant Complete with Controller, Battery Charger & Power Cord, TSMC, Rail Clamp & Cable Guide



Accessories & Spares	
Part Number	Description
BP37-003	Battery Charger
BP37-004	Power Cord UK
BP37-005	Power Cord EU
BP37-006	Power Cord USA
BP37-007	Power Cord AU
BP37-008	BP37 Thermal Support Matrix (Pack of 10) (<i>Matrix</i>)
BP37-009	BP37 Thermal Support Matrix Cable & Clamp (<i>TSMC</i>)
BP37-010	Fitting Kit for AHT Infant Transport Devices
BP37-014	Rail Clamp (LS Variant)
BP37-015	Cable Guide (LS Variant)
BP37-018	Universal Clamp



Advanced Healthcare Technology Ltd.

Theobalds Enterprise Centre, Platinum Way, Waltham Cross, Hertfordshire, EN8 8YD, United Kingdom

Distributed by:

*TransWarmer® is a Registered Trademark of Cooper Surgical inc.
Formula 1 is a Registered Trademark of Formula One Licencing B.V.*

*Baby Pod® are Registered Trademarks of AHT Ltd.
BabyPod 2™ & BP37™ are Trademarks of AHT Ltd.
BabyPod 2 is a patented product.*