

BP37™

Fitting Instructions

H - Variant



IMPORTANT:
TO USE THE BP37,
REFER TO APPROPRIATE IFU'S.




Parts Supplied:

Rail Clamp (Part No. BP37-014)


2 x M5 x12mm screws

Tools Required (not supplied):

4mm Allen Key

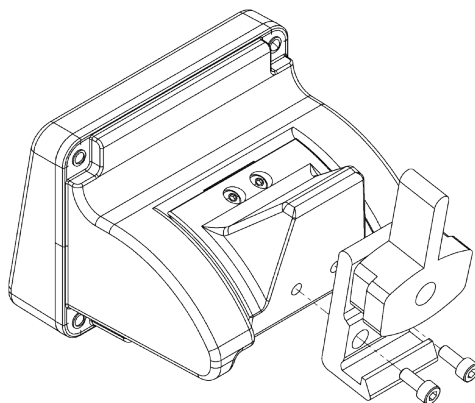

Attach Rail Clamp to Controller with supplied screws.

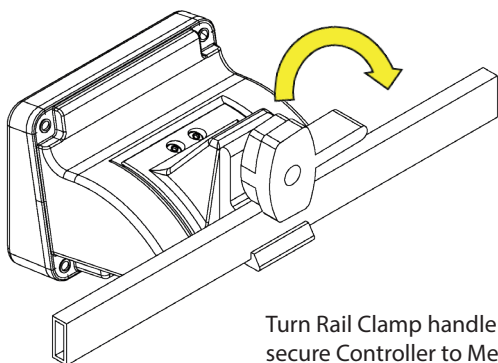

Turn Rail Clamp handle 90° to
secure Controller to Medi-Rail.

Reference: Use TSMC
(Part No. BP37-009-C)
with this controller.



Manufactured in France
on behalf of:



Advanced Healthcare Technology Ltd.
Theobalds Enterprise Centre,
Platinum Way, Waltham Cross,
Herts., EN8 8YD, United Kingdom
Tel: +44(0)1992 535933
Email: sales@babypod.com
www.babypod.com



Advanced Healthcare Technologies s.r.o.
Varšavská 715/36, 120 00 Prague 2, Czech Republic



BP37™ is a patented product.
BP37™ is a trademark of Advanced Healthcare Technology Ltd.